

MEDIKRO pro

PC-based Laboratory Spirometer

- ▶ All-in-one Laboratory Spirometer
- ▶ Embedded Ambient Condition Measurement
- ▶ Networking & Connectivity
- ▶ 3D Child Incentive



Product Code: M9488

Medikro® Pro represents the state of the art performance in PC-based spirometry. It is an accurate and precise diagnostic spirometer with embedded ambient condition measuring sensors for automatic real-time BTPS correction.

Pro utilizes the comprehensible Medikro® Spirometry Software (MSS). The system is a perfect fit for doctors, who require a precise spirometer with highly customizable software, clear visualization of results and excellent networking capabilities. With disposable Medikro® SpiroSafe flow transducers there is no need for bacterial filters or sterilization. This increases patient throughput and saves operational costs.



Pro Highlights



Freddy the Frog incentive

Medikro's Freddy the Frog is the only 3D incentive animation in the world. Freddy captures the attention of the children and encourages maximum performance. This improves the quality of spirometry results.

Cost-efficient

The combination of hardware and software features and the use of SpiroSafe disposables increases patient throughput and makes the system cost-efficient.



Patient safety

Disposable SpiroSafe flow transducers provide security in infection control. No exhaled air gets into the pressure tube or the device. This prevents patient-to-patient cross-contamination and protects the spirometer.



VCT technology

Vibration Control Tubing (VCT) cancels the signal noise caused by patient and tube movement. This stabilizes the signal and improves overall spirometry quality.



Automated realtime BTPS correction


The temperature, pressure and humidity sensors track environmental conditions automatically and enable true realtime BTPS correction.

Optional Bronchial Provocation

Use readymade bronchial provocation protocols, create them from the existing protocols or generate completely new ones.



Consumables & Accessories

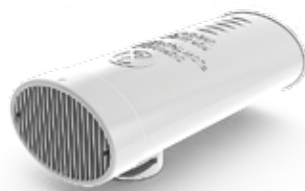

MEDIKRO® SpiroSafe

Product code: M9256



MEDIKRO® Calibration Syringe

Product code: M9474 (3L) / M9477 (1L)



🏠 Kellolahdentie 27, FI-70460 Kuopio, Finland
✉ P.O. Box 54, FI-70101, Kuopio, Finland
☎ +358 44 777 7002
📞 +358 17 283 3300
@ sales@medikro.com
🌐 www.medikro.com

Your local Medikro representative: