

# Medikro Pro

PRODUCT CODE  
M915

All-in-one PC-based laboratory spirometer

**Embedded Ambient Condition Measurement**

**Safe, disposable accessories**

**Networking & Connectivity**

**Interactive Child Incentive**



**MEDIKRO**  
Spirometry simplified



Medikro Pro represents the industry standard in PC-based spirometry. Highly accurate and precise diagnostic tool equipped with embedded sensors that measure ambient conditions. These sensors enable automatic real-time BTPS correction, which further enhances the accuracy of the results.

To conduct spirometry tests, Pro uses SpiroSafe flow transducers that are easy to handle and ensure protection of patients from cross-contamination. For added safety, Bacterial & Viral Filters are employed to protect operators. Both products are designed for single use, which makes it easy to accommodate more patients.

## PURCHASE BUNDLES

### Basic package - M985

Starter consumable pack

- M915 Medikro Pro Spirometer
- M9256-100 SpiroSafe Disposable Flow Transducer, 100 pieces
- M9264-200 Spare Pressure Tube with connectors for Flow Transducer, length 200 cm
- Nose Clip, 100 pieces

### Full package - M975

All components that the user needs

- M915 Medikro Pro Spirometer
- M9474 Calibration Syringe
- M9256-100 Disposable Flow Transducer, 100 pieces
- M9264-200 Spare Pressure Tube with connectors for Flow Transducer, length 200 cm
- MRN-01/100 Nose Clip, 100 pieces

## SERVICE AGREEMENT

Opt for a service agreement with Medikro to enhance your equipment's longevity and reliability.


### Indefinite Warranty

Service agreement customers receive a warranty extension as long as they continue annual maintenance, ensuring long-term value.

### Hassle-Free Maintenance

We handle scheduling and reminders, allowing you to focus on core responsibilities while maintaining optimal device performance.

## HIGHLIGHTS

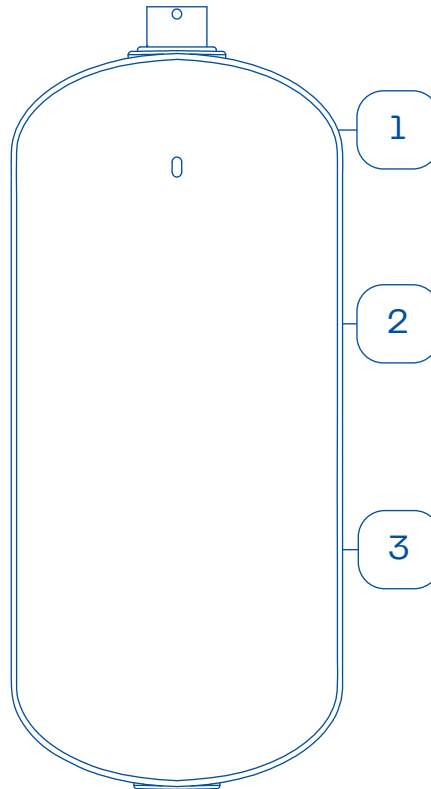

### VCT technology

Vibration Control Tubing (VCT) cancels the signal noise caused by patient and tube movement. This stabilizes the signal and improves overall spirometry quality.

### Automated realtime BTPS correction

The temperature, pressure and humidity sensors track environmental conditions automatically and enable true realtime BTPS correction.

### Optional Bronchial Provocation

Medikro Pro provides flexibility in bronchial provocation testing by allowing users to choose from ready-made protocols, create new ones based on existing templates, or generate completely novel protocols.

## CUSTOMERS' FEEDBACK

“ Medikro systems are reliable and affordable, their measurements precise and the software well thought out! ”

ENRIQUE CIMAS. GP FAMILY MEDICINE



Medikro Oy  
Pioneerinkatu 3  
70800 Kuopio, Finland

+358 44 777 7002  
www.medikro.com